diveMove: dive analysis in R

Sebastián P. Luque*

Contents

1	Introduction	1
2	Starting up	2
3	Reading Input Files	2
4	Extraction and Display of Information from TDR and TDRspeed Objects	3
5	ZOC and Wet/Dry period detection	4
6	Access to Elements from TDRcalibrate Objects	5
7	Speed Calibration	7
8	TDR dive and postdive statistics	8
9	Miscellaneous functions	10
10	Acknowledgements	10

1 Introduction

Dive analysis usually involves handling of large amounts of data, as new instruments allow for frequent sampling of variables over long periods of time. The aim of this package is to make this process more efficient for summarizing and extracting information gathered by time-depth recorders (TDRs, hereafter). The principal motivation for developing diveMove was to provide more flexibility during the various stages of

^{*}Contact: spluque@gmail.com. Comments for improvement are very welcome!

Table 1. Sample 1DA me structure.					
date	time	depth	light	temperature	speed
16/02/2004	14:30:00	12	200	8.4	1.44
16/02/2004	14:30:05	15	180	8.0	1.75
16/02/2004	14:30:10	19	170	7.6	1.99

Table 1. Sample TDR file structure.

analysis than that offered by popular commercial software. This is achieved by making the results from intermediate calculations easily accessible, allowing the user to make his/her own summaries beyond the default choices the package provides. The following sections of this vignette illustrate a typical work flow during analysis of TDR data, using the sealMK8 data available in diveMove as an example.

2 Starting up

As with other packages in R, to use the package we load it with the funtion library:

> library(diveMove)

This makes the objects in the package available in the current R session. A short overview of the most important functions can be seen by running the examples in the package's help page:

> example(diveMove)

3 Reading Input Files

Input files must be simple, comma-delimited text files¹. The order of columns is not significant, as the column numbers indicating the variables of interest can be supplied as arguments. Table 1 shows the file structure that readTDR assumes by default, which is a standard structure of files from common TDR models.

Depending on the TDR model, speed may be omitted. Currently, light, temperature and any other variables beyond column 6 are ignored.

To read the file into R, use the function readTDR:

¹The extension does not matter, but conventionally these files have a .csv extension

4 Extraction and Display of Information from TDR and TDRspeed Objects

Read the help page for readTDR using ?readTDR following common R help facilities. Thus, data could have been subsampled at a larger interval than that in the original file, so that the time interval between readings is 10 s:

But since the original 5 s interval (which is the default value for *subsamp*) is what will be used for the subsequent sections, it is recreated it here:

The format in which date and time should be interpreted can be controlled with the argument *dtformat*. If the data are already available in the R session, e.g. as a data frame, then the function createTDR can be used to convert it to a form that facilitates further analyses.

Both of these functions store the data in an object of class TDR or TDRspeed, which hold information on the source file and sampling interval, in addition to the variables described above. Which of these objects is created is determined by the name of the input file. All files should contain the letter sequence "mk" followed by a number, as these correspond to the names of common TDR models. If the number following this sequence is 8, then a TDRspeed object is created, otherwise the function returns a TDR object.

4 Extraction and Display of Information from TDR and TDRspeed Objects

For convenience, extractor methods are available to access the different slots from objects of these classes. The standard *show* method will display the usual overview information on the object:

> sealX

Time-Depth Recorder data -- Class TDRspeed object

Source File : sealMK8.csv

Sampling Interval (s) : 5 Number of Samples : 34199

Total Duration (d) : 1.979

Other extractor methods are named after the component they extract: getTime, get-Depth, getSpeed, and getDtime, where the latter extracts the sampling interval. The plot method brings up a plot of the data covering the entire record, although a tcltk widget provides controls for zooming and panning to any particular time window. Alernatively, the underlying function plotDive provides the same functionality, but takes separate time and depth arguments, rather than a TDR object.

At any time, TDR objects can be coerced to a simple data frame, which can later be exported or manipulated any other way:

```
> sealX.df <- as.data.frame(sealX)
> head(sealX.df)
[1] 5
```

5 Zero-Offset Depth Correction and Summary of Wet/Dry Periods

One the first steps of dive analysis involves correcting depth for shifts in the pressure transducer, so that surface readings correspond to the value zero. Although some complex algorithms exist for detecting where these shifts occur in the record, the shifts remain difficult to detect and dives are often missed, which a visual examination of the data would have exposed. The trade off is that visually zero-adjusting depth is tedious, but the advantages of this approach far outweigh this cost, as much insight is gained by visually exploring the data. Not to mention the fact that obvious problems in the records are more effectively dealt with in this manner.

That personal rant aside, zero offset correction (ZOC) is done in diveMove using the function zoc. However, a more efficient method of doing this is by using the calibrat-eDepth function, which takes a TDR object (or inheriting from it) to perform three basic tasks. The first is to ZOC the data, using the tcltk package to be able to do it interactively:

> dcalib <- calibrateDepth(sealX)</pre>

This command brings up a plot with tcltk controls allowing to pan and zoom in or out of the data, as well as adjustment of the depth scale. Thus, an appropriate time window with a unique surface depth value can be displayed. It is important to make the display such that the depth scale is small enough to allow the resolution of the surface value with the mouse. Clicking on the ZOC button waits for two clicks:

1. the coordinates of the first click define the starting time for the window to be ZOC'ed, and the depth corresponding to the surface,

2. the second click defines the end time for the window (only the x coordinate has any meaning).

This procedure can be repeated as many times as needed. If there is any overlap between time windows, then the last one prevails. However, if the offset is known a priori, there is no need to go through all this procedure, and the value can be provided as the argument *offset* to calibrateDepth.

Once depth has been ZOC'ed, calibrateDepth will identify dry and wet periods in the record. Wet periods are those where a depth reading is available, dry periods are those without a depth reading. Records often have abherrant missing depth that should not be considered dry periods, as they are often of very short duration. Likewise, there may be periods of wet activity that are too short to be compared with other wet periods. This can be controlled by setting the arguments landerr and seaerr.

Finally, calibrateDepth identifies all dives in the record, according to a minimum depth criteria given as its *divethres* argument. The result (value) of this function is an object of class *TDRcalibrate*, where all the information obtained during the tasks described above are stored. Again, an appropriate *show* method is available to display a short overview of such objects:

> dcalib

Depth calibration -- Class TDRcalibrate object

Source file : sealMK8.csv

Number of dry phases : 4

Number of aquatic phases : 3

Number of dives detected : 317

Dry threshold used (s) : 70

Aquatic theshold used (s) : 3610

Dive threshold used (s) : 4

Speed calibration coefficients : a = 0; b = 1

6 Access to Elements from TDRcalibrate Objects

Extractor methods are also available to access the information stored in *TDRcalibrate* objects. These include: getTDR, getGAct, getDAct, getDPhaseLab, and getSpeedCoefs. These are all generic functions² that access the (depth) calibrated TDR object, details from wet/dry periods, dives, dive phases, and speed calibration coefficients (see Section 7), respectively. Below is a short explanation of these methods.

²A few of them with more than one method

getTDR This method simply takes the TDR calibrate object as its single argument and extracts the TDR object³:

> getTDR(dcalib)

Time-Depth Recorder data -- Class TDRspeed object

Source File : sealMK8.csv

Sampling Interval (s) : 5 Number of Samples : 34199

Sampling Begins : 2002-01-05 11:32:00 Sampling Ends : 2002-01-07 11:01:50

Total Duration (d) : 1.979

getGAct There are two methods for this generic, allowing access to a list with details about all wet/dry periods found. One of these extracts the entire *list* (output omitted for brevity):

> getGAct(dcalib)

The other provides access to particular elements of the *list*, by their name. For example, if we are interested in extracting only the vector that tells us to which period number every row in the record belongs to, we would issue the command:

> getGAct(dcalib, "phase.id")

Other elements that can be extracted are named "trip.act", "trip.beg", and "trip.end", and can be extracted in a similar fashion. These elements correspond to the activity performed for each reading (see ?detPhase for a description of the labels for each activity), the beginning and ending time for each period, respectively.

getDAct This generic also has two methods; one to extract an entire data frame with details about all dive and postdive periods found (output omitted):

> getDAct(dcalib)

The other method provides access to the columns of this data frame, which are named "dive.id", "dive.activity", and "postdive.id". Thus, providing any one of these strings to getDAct, as a second argument will extract the corresponding column.

getDPhaseLab This generic function extracts a factor identifying each row of the record to a particular dive phase (see ?detDive for a description of the labels of the factor identifying each dive phase). Two methods are available; one to extract the entire factor, and the other to select particular dive(s), by its (their) number, respectively (output omitted):

³In fact, a *TDRspeed* object in this example

```
> getDPhaseLab(dcalib)
> getDPhaseLab(dcalib, 20)
> dphases <- getDPhaseLab(dcalib, c(100:300))</pre>
```

The latter method is useful for visually inspecting the assignment of points to particular dive phases. Before doing that though, this is a good time to introduce another generic function that allows the subsetting of the original TDR object to a single a dive or group of dives' data:

As can be seen, the function takes a TDRcalibrate object and a vector indicating the dive numbers to extract, and returns a TDR object containing the subsetted data. Once a subset of data has been selected, it is possible to plot them and pass the factor labelling dive phases as the argument phaseCol to the plot method⁴:

```
> plot(subSealX, phaseCol = dphases)
```

7 Speed Calibration

Calibration of speed readings is done using the principles described in Blackwell (1999) and Hindell et al. (1999). The function calibrateSpeed performs this operation⁵, and allows the selection of the particular subset of the data that should be used for the calibration:

```
> vcalib <- calibrateSpeed(dcalib, calType = "pooled")
> vcalib
```

⁴The function that the method uses is actually plotDive, so all the possible arguments can be studied by reading the help page for plotDive

⁵CAUTION: This implementation is experimental, and may give unexpected results.

Depth calibration -- Class TDRcalibrate object

Source file : sealMK8.csv

Number of dry phases : 4

Number of aquatic phases : 3

Number of dives detected : 317

Dry threshold used (s) : 70

Aquatic theshold used (s) : 3610

Dive threshold used (s) : 4

Speed calibration coefficients : a = 0.4; b = 0.64

A side effect of such a call is the production of a plot displaying the quantile regression fit for the three phases (Figure 1). This can be displayed on the current device, or sent to a postscript file, using postscript=TRUE in the call, for a higher quality representation.

The default (calType="pooled") is to use data from the descent and ascent phases of all dives, but possible values also include "descent", "ascent", and "none". Because the function produces three plots of speed vs. rate of depth change, the latter is useful in cases where speed does not need any calibration, but inspection of the plots is desired. Finer control is possible by the use of arguments type, which controls whether descent or ascent readings that are shared with the bottom phase of the dive should be included or not, and bad, which controls minimum speeds and rates of depth change through which the calibration line should be drawn. Finally, a maximum depth threshold can be supplied as the argument z, so that only data from dives where maximum depth was greater than this value are included in the construction of the calibration line.

If the calibration coefficients from the implicit quantile regression are known a priori, then these can be supplied to the function via its *coefs* argument. In this case, no plots are created.

8 TDR dive and postdive statistics

Once data have been calibrated and the record broken up at "trip" and "dive" scales, obtaining dive statistics is a trivial call to function diveStats:

```
> dives <- diveStats(vcalib)
> head(dives, 3)
```

```
begdesc enddesc

1 2002-01-05 12:20:10 2002-01-05 12:20:15

2 2002-01-05 21:19:40 2002-01-05 21:20:20

3 2002-01-05 21:22:10 2002-01-05 21:23:15
```

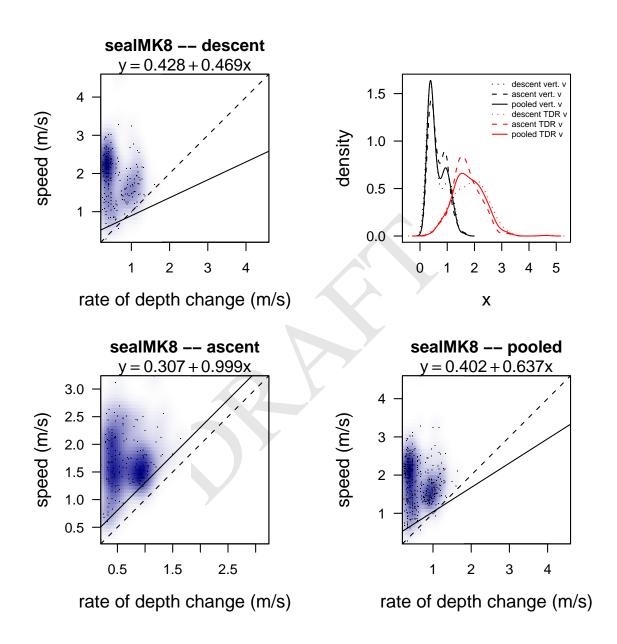


Figure 1. Example speed calibration lines, dividing dives into descent, ascent, or pooling both phases from a TDR record.

		begasc	desctim	botttim	asctim o	descdist
1	2002-01-05	12:20:20	7.5	5	7.5	6
2	2002-01-05	21:20:40	42.5	20	47.5	26
3	2002-01-05	21:23:50	67.5	35	72.5	63
	bottdist as	scdist de	sc.tdist	desc.mea	an.speed	desc.angle
1	0	6	22.44		4.488	15.51
2	3	29	100.07		2.502	15.06
3	8	67	107.84		1.659	35.75
	bott.tdist	bott.mean	n.speed a	asc.tdist	asc.me	an.speed
1	15.22		3.043	18.04	<u> </u>	3.609
2	53.96		2.698	71.78	3	1.595
3	56.11		1.603	98.09)	1.401
	asc.angle	divetim m	axdep pos	stdive.du	ır postd:	ive.tdist
1	19.42	20	6	3234	! 5	50445.70
2	23.83	110	29	3	35	16.85
3	43.08	175	67	7	' 5	58.18
	postdive.me	ean.speed				
1		1.5652				
2		0.4815				7
3		0.7758				

The function takes a single argument: an object of class *TDRcalibrate*, and returns a data frame with one row per dive in the record, with a suite of basic dive statistics in each column. Please consult ?diveStats for an explanation of each of the variables estimated, although the names of the output data frame should be self explanatory. These variables are thus available for calculating any other derived values, by extracting them using the standard R subscripting facilities.

9 Miscellaneous functions

Other functions are included for handling location data, and these are readLocs, aust-Filter, and distSpeed. These are useful for reading, filtering, and summarizing travel information. For extensive animal movement analyses, refer to package timeTrack available at http://staff.acecrc.org.au/~mdsumner/Rutas/.

10 Acknowledgements

Invaluable input and help during development of this package has been offered by John P.Y. Arnould, and regular contributors to R-help.

References

Blackwell, S. (1999). A method for calibrating swim-speed recorders. *Marine Mammal Science*, 15(3):894–905.

Hindell, M., McConnell, B., Fedak, M., Slip, D., Burton, H., Reijnders, P., and McMahon, C. (1999). Environmental and physiological determinants of successful foraging by naive southern elephant seal pups during their first trip to sea. *Canadian Journal of Zoology*, 77:1807–1821.



diveMove

November 15, 2006

R topics documented:

	ustFilter	1
	calibrateDepth	3
	letDive	4
	letPhase	5
	listSpeed	7
	liveMove-internal	8
	liveMove-package	9
	liveStats	10
	eadLocs	12
	eadTDR	13
	qPlot	14
	ealMK8	15
	loSpeedCalib	16
	TDRcalibrate-class	17
	TDR-class	20
	toc	22
Index		24

austFilter

Filter satellite locations

Description

Apply a three stage algorithm to eliminate erroneous locations, based on the procedure outlined in Austin et al. (2003).

Usage

```
austFilter(time, lon, lat, id=gl(1, 1, length(time)),
   speedthres, distthres, window=5)
grpSpeedFilter(x, speedthres, window=5)
rmsDistFilter(x, speedthres, window=5, distthres)
```

2 austFilter

Arguments

time	POSIXct object with dates and times for each point.
lon	numeric vectors of longitudes, in decimal degrees.
lat	numeric vector of latitudes, in decimal degrees.
id	a factor grouping points in different categories (e.g. individuals).
speedthres	speed threshold (m/s) above which filter tests should fail any given point.
distthres	distance threshold above which the last filter test should fail any given point.
window	integer indicating the size of the moving window over which tests should be carried out.
х	3-column matrix with column 1: POSIXct vector; column 2: numeric longitude vector; column 3: numeric latitude vector.

Details

These functions implement the location filtering procedure outlined in Austin et al. (2003). grpSpeedFilter and rmsDistFilter can be used to perform only the first stage or the second and third stages of the algorithm on their own, respectively. Alternatively, the three filters can be run sequentially using austFilter.

The first stage of the filter is an iterative process which tests every point, except the first and last two, for travel velocity relative to the preceeding/following two points. If all these four speeds are greater than the specified threshold, the point is marked as failing the first stage. In this case, the next point is tested, removing the failing point from the set of test points.

The second stage runs McConnell et al. (1992) algorithm, which tests all the points that passed the first stage, in the same manner as above. The root mean square of all four speeds is calculated, and if it is greater than the specified threshold, the point is marked as falling the second stage.

The third stage is run simultaneously with the second stage, but if the mean distance of all four pairs of points is greater than the specified threshold, then the point is marked as failing the third stage.

Value

A matrix with three columns of logical vectors with values TRUE for points that failed each stage. Results from each filter are presented independently of the others; i.e. points marked as failing one filter are not necessarily marked as failing the next one.

Warning

This function applies McConnell et al.'s filter as described in Austin et al. (2003), but other authors may have used the filter differently. Austin et al. (2003) have apparently applied the filter in a vectorized manner. It is not clear from the original paper whether the filter is applied iteratively or in a vectorized fashion, so authors may be using it inconsistently.

Author(s)

Sebastian P. Luque (spluque@gmail.com) and Andy Liaw.

calibrateDepth 3

References

McConnell BJ, Chambers C, Fedak MA. 1992. Foraging ecology of southern elephant seals in relation to bathymetry and productivity of the Southern Ocean. *Antarctic Science* 4:393-398.

Austin D, McMillan JI, Bowen D. 2003. A three-stage algorithm for filtering erroneous Argos satellite locations. *Marine Mammal Science* 19: 371-383.

See Also

```
distSpeed
```

calibrateDepth Calibrate and build a "TDRcalibrate" object

Description

These functions create a "TDRcalibrate" object which is necessary to obtain dive summary statistics.

Usage

Arguments

```
x an object of class TDR for calibrateDepth, and an object of class TDR calibrate-class for calibrateSpeed.

landerr, seaerr arguments to detPhase.

divethres argument to zoc.

type, calType, bad, z, filename further arguments for .getSpeedCalib and doSpeedCalib.

coefs known speed calibration coefficients from quantile regression as a vector of length 2 (intercept, slope). If provided, these coefficients are used for calibrating speed, ignoring all other arguments, except x.

... argument passed to doSpeedCalib.
```

Details

These functions are really wrappers around functions that are usually called in sequence, so they provided an abbreviated method for running them together during analyses. See the functions in the 'See Also' section for more details.

calibrateDepth performs zero-offset correction of depth, wet/dry phase detection, and detection of dives, as well as proper labelling of the latter.

calibrateSpeed calibrates speed readings.

detDive. 4

Value

```
An object of class TDRcalibrate-class
```

Author(s)

Sebastian P. Luque (spluque@gmail.com)

See Also

detPhase, detDive, doSpeedCalib, zoc, for the underlying functions.

detDive

Detect dives from depth readings

Description

Identify dives in TDR records based on a dive threshold.

Usage

```
detDive(zdepth, act, divethres=4, ...)
labDive(act, string, interval)
labDivePhase(x, diveID)
```

Arguments

vector of zero-offset corrected depths. zdepth

factor as long as depth coding activity, with levels specified as in detPhase. divethres threshold depth below which an underwater phase should be considered a dive. string a character belonging to a level of act to search for and label sequentially.

interval,

the sampling interval in seconds.

a class 'TDR' object

numeric vector indexing each dive (non-dives should be 0) diveID

Details

emph { detDive} detects a dive whenever the zero-offset corrected depth in an underwater phase is below the supplied dive threshold. The adjustment is done only for phases of at-sea activity, completely ignoring phases with other activity.

emph{labDive} assigns a unique number to each dive along a vector of depths, and equally numbering the subsequent postdive interval.

emph{labDivePhase} labels each row identifying it with a portion of the dive.

detPhase 5

Value

A data frame with the following elements for detDive

```
dive.id numeric vector numbering each dive in the record.
```

factor with levels 'L', 'W', 'U', 'D', and 'Z', see $\mathtt{detPhase}$. All levels may be represented.

postdive.id numeric vector numbering each postdive interval with the same value as the preceding dive.

labDive returns a matrix with as many rows as its first two arguments with two columns: dive.id, and postdive.id, each one sequentially numbering each dive and postdive period.

labDivePhase returns a factor with levels "D", "DB", "B", "B", "A", "DA", and "X", breaking the input into descent, descent/bottom, bottom/ascent, ascent, and non-dive, respectively.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

See Also

```
detPhase, zoc
```

Detect phases of activity from depth readings

Description

detPhase

Functions to identify sections of a TDR record displaying one of three possible activities: on-land, at-sea, and at-sea leisure.

Usage

```
detPhase(time, depth, landerr, seaerr, ...)
getAct(time, act, interval)
```

Arguments

landerr

time POSIXct object with date and time for all depths.

depth numeric vector with depth readings.

land error threshold in seconds. On-land phases shorter than this threshold will

be considered as at-sea.

seaerr at-sea leisure threshold in seconds. At-sea phases shorter than this threshold will

be considered as at-sea leisure.

act A numeric vector indicating the activity for every element of time.

interval, .

sampling interval in seconds.

6 detPhase

Details

detPhase first creates a factor with value 'L' (on-land) for rows with NAs for depth and value 'W' (at-sea) otherwise. It subsequently calculates the duration of each of these phases of activity. If the duration of an on-land phase ('L') is less than landerr, then the values in the factor for that phase are changed to 'W' (at-sea). The duration of phases is then recalculated, and if the duration of a phase of at-sea activity is less than seaerr, then the corresponding value for the factor is changed to 'Z' (at-sea leisure). The durations of all phases are recalculated a third time to provide final phase durations.

getAct takes a factor indicating different activity phases, their associated time, and the sampling interval to return a factor uniquely identifying each phase of activity, i.e. labelling them. In addition, it returns the duration of each phase, and their beginning and end times.

Value

A list with components; the first 4 are returned by detPhase and the rest by getAct:

phase.id	numeric vector identifying each activity phase, starting from $\boldsymbol{1}$ for every input record.
trip.act	factor with levels 'L' indicating land, 'W' indicating at-sea, 'U' for underwater (above dive criterion), 'D' for diving, 'Z' for at-sea leisure animal activities. Only 'L', 'W', and 'Z' are actually represented.
trip.beg	a POSIXct object as long as the number of unique activity phases identified, indicating the start times for each activity phase.
trip.end	a POSIXct object as long as the number of unique activity phases identified, indicating the end times for each activity phase.
time.br	a factor dividing the factor act in phases.
time.peract	duration of each phase defined by time.br.
beg.time	POSIXct object; beginning time for each phase.
end.time	POSIXct object; ending time for each phase.

Author(s)

Sebastian P. Luque (spluque@gmail.com) and Andy Liaw.

See Also

detDive

distSpeed 7

distSpeed Calculate distance and speed between locations
distSpeed Calculate distance and speed between locations

Description

Calculate distance, time difference, and speed between pairs of points defined by latitude and longitude, given the time at which all points were measured.

Usage

```
distSpeed(pt1, pt2, speed=TRUE)
track(txy, id=q1(1, nrow(txy)), subset)
```

Arguments

pt1	a matrix or data frame with three columns; the first a POSIXct object with dates and times for all points, the second and third numeric vectors of longitude and latitude for all points, respectively, in decimal degrees.
pt2	a matrix with the same structure as pt1.
speed	logical; should speed between points be calculated?
txy	a data frame with a POSIXct object in its first column, lon and lat in second and third column, respectively.
id	a factor dividing the data in txy into distinct groups.
subset	a logical expression indicating the rows to be analyzed, in terms of elements of txy .

Details

pt1 and pt2 may contain any number of rows. track is essentially a wrapper for distSpeed, taking a data frame, assumed to be ordered chronologically, and calculations are done between all successive rows.

Value

For distSpeed, a matrix with three columns: distance (km), time difference (h), and speed (m/s). For track, a data frame with an id column and the same columns as in distSpeed.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

8 diveMove-internal

diveMove-internal Internal diveMove Functions

Description

Functions used for very particular tasks within larger functions in diveMove

Usage

```
cutDive(x)
.descAsc(x, phase, type=c("all", "strict"), interval, z=0)
.diveIndices(diveID, diveNo)
.getInterval(time)
.getSpeedCalib(time, zdepth, speed, dives, phase, ...)
.getSpeedStats(x, vdist)
```

Arguments

х	a single dive's data; for <code>.cutDive</code> : a 2-col matrix with subscript in original TDR object and non NA depths. For <code>.descAsc</code> : a 4-col matrix with dive id, time, depth, and speed. For <code>.getSpeedStats</code> : a 3-col matrix with time, depth, and speed.
time	POSIXct object representing time.
dtformat	A string to interpret date and time (see strptime.
phase	factor labelling each row for its phase in dive.
type	string indicating whether all points belonging to descent/ascent should be included ("all"), or points shared with bottom phase should be excluded ("srict").
interval	sampling interval in seconds.
z	minimum depth differences to use.
zdepth	zero-offset corrected depth m.
speed	speed in m/s. For doSpeedCalib: uncalibrated speeds; ignored if ${\tt calType}$ is "none".
dives	3-col data.frame with dive id (numeric), activity (factor), and postdive id (numeric).
	arguments to pass to .descAsc (type, interval, and z).
vdist	vertical distance travelled during ascent or descent.
diveID	Numeric vector of all dive and non dive IDs.
diveNo	Numeric vector of unique dive IDs to index in ${\tt diveID}$.

diveMove-package 9

Details

These functions are not meant to be called directly by the user, as he/she could not care less (right?). This may change in the future.

.getSpeedCalib extracts the rates of descent and ascent with associated mean speed during descent and ascent phases, respectively and returns a list that is later manipulated by doSpeedCalib to calibrate speed. The speed used for each rate of depth change corresponds to the speed read for the last point, assuming that each speed reading is the average speed for the last measurement interval.

Value

.getSpeedCalib: A list with two elements (named "descent" and "ascent"). Each element is a 2-column matrix with rate of depth change in the first column, and speed in the second, corresponding to the descent phase of each dive.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

diveMove-package Time depth recorder analysis

Description

This package is a collection of functions for visualizing, and analyzing depth and speed data from time-depth recorders *TDRs*. These can be used to zero-offset correct depth, calibrate speed, and divide the record into different phases, or time budget. Functions are provided for calculating summary dive statisfics for the whole record, or at smaller scales within dives.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

See Also

A vignette with a guide to this package is available by doing 'vignette("diveMove")'. TDR-class, calibrateDepth, calibrateSpeed, attendance, stampDive.

Examples

10 diveStats

```
## Not run: dcalib <- calibrateDepth(sealX) # interactively
(dcalib <- calibrateDepth(sealX, offset=3)) # zero-offset correct at 3 m
## Not run:
## plot all readings and label them with the phase of the record
## they belong to, excluding surface readings
plot(dcalib, labels="phase.id", surface=FALSE)
## plot the first 300 dives, showing dive phases and surface readings
plot(dcalib, diveNo=seg(300), labels="dive.phase", surface=TRUE)
## End(Not run)
## calibrate speed (using default remaining arguments)
(vcalib <- calibrateSpeed(dcalib))
## Obtain dive statistics for all dives detected
dives <- diveStats(vcalib)
head (dives)
## Attendance table
att <- attendance(vcalib, FALSE) # taking trivial aquatic activities into account
att <- attendance (vcalib, TRUE) # ignoring them
## Add trip stamps to each dive
stamps <- stampDive(vcalib)
sumtab <- data.frame(stamps, dives)
head(sumtab)
```

diveStats

Per-dive statistics

Description

Calculate dive statistics in TDR records

Usage

```
diveStats(x)
getDive(x, interval, speed=FALSE)
stampDive(x, ignoreZ=TRUE)
```

Arguments

interval

a TDR calibrate-class object for diveStats and stampDive, a data v frame containing a single dive's data.

sampling interval for interpreting x.

logical; should speed statistics be calculated? speed

logical indicating whether trips should be numbered considering all aquatic acignoreZ

tivities ("W" and "Z") or ignoring "Z" activities.

diveStats 11

Details

diveStats calculates various dive statistics based on time and depth for an entire TDR record. getDive obtains these statistics from a single dive, and stampDive stamps each dive with associated trip information.

Value

A data.frame with one row per dive detected (durations are in s, and linear variables in m):

begdesc	A POSIXct object, specifying the start time of each dive.		
enddesc	A POSIXct object, as begdesc indicating descent's end time.		
begasc	A POSIXct object, as begdesc indicating the time ascent began.		
desctim	descent duration of each dive.		
botttim	bottom duration of each dive.		
asctim	ascent duration of each dive.		
descdist	numeric vector with descent depth.		
bottdist	numeric vector with the sum of absolute depth differences while at the bottom of each dive; measure of amount of "wiggling" while at bottom.		
ascdist	numeric vector with ascent depth.		
desc.tdist	numeric vector with descent total distance, estimated from speed.		
desc.mean.sp	eed		
	numeric vector with descent mean speed.		
desc.angle	numeric vector with descent angle.		
bott.tdist	numeric vector with bottom total distance, estimated from speed.		
bott.mean.sp	bott.mean.speed		
	numeric vector with bottom mean speed.		
asc.tdist	numeric vector with ascent total distance, estimated from speed.		
asc.mean.spe	ed		
	numeric vector with ascent mean speed.		
asc.angle	numeric vector with ascent angle.		
divetim	dive duration.		
maxdep	numeric vector with maximum depth.		
postdive.dur	postdive duration.		
postdive.tdist			
	numeric vector with postdive total distance, estimated from speed.		
postdive.mea	n . speed numeric vector with postdive mean speed.		

The number of columns depends on the value of speed.

stampDive returns a data frame with trip number, trip type, and start and end times for each dive.

12 readLocs

Author(s)

Sebastian P. Luque (spluque@gmail.com)

See Also

```
detPhase, zoc, TDRcalibrate-class
```

readLocs Read comma-delimited file with location data

Description

Read a comma delimited (*.csv) file with (at least) time, latitude, longitude readings.

Usage

```
readLocs(file, loc.idCol, idCol, dateCol, timeCol=NULL,
    dtformat="%m/%d/%Y %H:%M:%s", tz="GMT",
    classCol, lonCol, latCol, alt.lonCol=NULL, alt.latCol=NULL)
```

Arguments

file	A string indicating the name of the file to read. Provide the entire path if the file is not on the current directory.
loc.idCol	Column number containing location ID.
idCol	Column number containing an identifier for locations belonging to different groups.
dateCol	Column number containing dates, and, optionally, times.
timeCol	Column number containing times.
dtformat	A string, specifying the format in which the date and time columns, when pasted together, should be interpreted (see strptime) in file.
tz	A string indicating the time zone for the date and time readings.
lonCol	Column number containing longitude readings.
latCol	Column number containing latitude readings.
classCol	Column number containing the ARGOS rating for each location.
alt.lonCol	Column number containing alternative longitude readings.
alt.latCol	Column number containing alternative latitude readings.

Details

The file must have a header row identifying each field, and all rows must be complete (i.e. have the same number of fields). Field names need not follow any convention.

readTDR 13

Value

A data frame.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

readTDR

Read comma-delimited file with TDR data

Description

Read a comma delimited (*.csv) file containing time-depth recorder (TDR) data from various TDR models. Models supported are MK5, MK7, and MK8 Wildlife Computers instruments. Return a TDR or TDRspeed object. buildTDR creates an object of one of these classes from other objects in the session.

Usage

```
readTDR(file, dateCol=1, timeCol=2, depthCol=3, speedCol=6,
    subsamp=5, dtformat="%d/%m/%Y %H:%M:%S", tz="GMT")
    createTDR(time, depth, speed, dtime, file)
```

Arguments

file	A string indicating the path to the file to read.
dateCol	Column number containing dates, and optionally, times.
timeCol	Column number with times.
depthCol	Column number containing depth readings.
speedCol	Column number containing speed readings.
subsamp	Subsample rows in file with subsamp interval, in s.
dtformat	A string, specifying the format in which the date and time columns, when pasted together, should be interpreted (see ${\tt strptime}$).
tz	A string indicating the time zone assumed for the date and time readings.
time	a ${\tt POSIXct}$ object with date and time readings for locations.
depth	numeric vector with depth readings.
speed	optional numeric vector with speed readings.
dtime	sampling interval used in seconds.

14 rqPlot

Details

The file name must contain the adjacent letter "mk" somewhere to be able to identify the TDR model. If the number following these letters is an 8, then a column for speed readings is expected, in addition to depth.

The file must have a header row identifying each field, and all rows must be complete (i.e. have the same number of fields). Field names need not follow any convention. However, depth and speed should preferably be given in m, and $m \cdot s^{-1}$ for further analyses.

Value

An object of class 'TDR' or 'TDRspeed'.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

Examples

rqPlot

Plot of quantile regression for speed calibrations

Description

Plot of quantile regression for assessing quality of speed calibrations

Usage

```
rgPlot(rdepth, speed, rqFit, main="qtRegression",
    xlab="rate of depth change (m/s)", ylab="speed (m/s)",
    colramp=colorRampPalette(c("white", "darkblue")))
```

Arguments

speed in m/s.

rdepth numeric vector with rate of depth change.

 $\verb"rqFit" object of class "rq" representing a quantile regression fit of rate of depth change$

on mean speed.

main string; title prefix to include in ouput plot.

xlab, ylab axis labels.

colramp function taking an integer n as an argument and returning n colors.

sealMK8

Details

The dashed line in the plot represents a reference indicating a one to one relationship between speed and rate of depth change. The other line represent the quantile regression fit.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

See Also

```
doSpeedCalib, diveStats
```

sealMK8

Sample TDR data from a fur seal

Description

This data set is meant to show the organization a TDR *.csv file must have in order to be used as input for readTDR.

Format

A comma separated value (csv) file with 69560 TDR readings with the following columns:

date date

time time

depth depth in m

light light level

temperature temperature in C

speed speed in m/s

Details

The data is a subset of an entire TDR record, so it is not meant to make any inferences from this particular individual/deployment.

Author(s)

Sebastian P. Luque (spluque@gmail.com)

Source

Sebastian P. Luque, Christophe Guinet, John P.Y. Arnould

See Also

readTDR

16 doSpeedCalib

doSpeedCalib Calibration of TDR speed

Description

Calibrate speed readings from a TDR, based on the principles outlined in Blackwell et al. (1999).

Usage

Arguments

rates	two-element list corresponding to descent and ascent phases of dives, respectively. Each element should be a 3-column matrix with dive id, rate of depth
	change, and mean speed.
speed	numeric vector with uncalibrated speeds.
calType	string specifying the type of calibration to perform. It should be one of "descent", "ascent", or "pooled".
bad	vector of length 2 indicating values for rate of depth change and mean speed, respectively, below which data should be excluded to build the calibration curve.
filename	string indicating name of file to use as base name for the output postscript file.
postscript	logical; whether output plot to eps.
	arguments passed to rqPlot; currently, xlab, ylab, and colramp.

Details

Provide calibrated speeds in a TDR record, using the quantile regression of speed on rate of depth change, based on principles outlined in Blackwell et al. (1999). Choice of calibrating against pooled, or descent phases.

The function takes the rates of depth change and speed, for each phase of the dive separately or combined (based on the value of calType). It subsequently fits a quantile regression through the second percentile of the distribution of speed conditional on rate of depth change. The calibrated speed is $s_c = \frac{s_u - a}{b}$, where s_c is the calibrated speed, s_u is the uncalibrated speed, and a and b are the intercept and slope of the quantile regression, respectively.

Value

If calType is not "none", a list of two elements:

coefficients numeric vector of length two with the intercept and the slope of the quantile regression defining the calibration curve.

corrSpeed numeric vector as long as speed with the calibrated speeds.

A plot (possibly via postscript, depending on the value of postscript argument) of the calibration lines for all possible cases: "descent", "ascent", and "pooled", is created as a side effect.

TDRcalibrate-class 17

Author(s)

```
Sebastian P. Luque (spluque@gmail.com)
```

References

Blackwell, S. (1999) A method for calibrating swim-speed recorders. Marine Mammal Science 15(3): 894-905.

See Also

```
TDRcalibrate-class, rgPlot
```

TDRcalibrate-class Class "TDRcalibrate" for dive analysis

Description

This class holds information produced at various stages of dive analysis. Methods are provided for extracting data from each slot.

Usage

```
## S4 method for signature 'TDRcalibrate, logical':
attendance(obj, ignoreZ)
## S4 method for signature 'TDRcalibrate, numeric,
    missing':
extractDive(obj, diveNo)
## S4 method for signature 'TDRcalibrate, missing':
getDAct(x)
## S4 method for signature 'TDRcalibrate, character':
getDAct(x, y)
## S4 method for signature 'TDRcalibrate, missing':
getDPhaseLab(x)
## S4 method for signature 'TDRcalibrate, numeric':
getDPhaseLab(x, diveNo)
## S4 method for signature 'TDRcalibrate, missing':
getGAct(x)
## S4 method for signature 'TDRcalibrate, character':
getGAct(x, y)
## S4 method for signature 'TDRcalibrate, missing':
plot(x, diveNo=seg(unique(getDAct(x, "dive.id"))),
     labels="phase.id", surface=FALSE, ...)
```

18 TDRcalibrate-class

Arguments

x, obj	"TDRcalibrate" object.		
ignoreZ	logical indicating whether to ignore trivial aquatic periods.		
diveNo	numeric vector with dive numbers to plot.		
labels	one of "phase.id" or "dive.phase", specifying whether to label observations based on the gross phase ID of the "TDR" object, or based on each dive phase, respectively.		
surface	logical indicating whether to plot surface readings.		
	further arguments to plotDive.		
У	string; "dive.id", "dive.activity", or "postdive.id" in the case of getDAct, to extract the numeric dive ID, the factor identifying dive phases in each dive, or the numeric postdive ID, respectively. In the case of getGAct it should be one of "phase.id", "trip.act", "trip.beg", or "trip.nd", to extract the numeric phase ID for each observation, a factor indicating what major activity the observation corresponds to, or the beginning and end times of each phase in the record, respectively.		

Details

This is perhaps the most important class in diveMove, as it holds all the information necessary for calculating requested summaries for a TDR.

The tdr slot contains the time, zero-offset corrected depth, and possibly calibrated or uncalibrated speed. See readTDR and the accessor function getTDR for this slot. Convenient access to each vector in this slot is available through getTime, getDepth, and getSpeed.

The slot gross.activity holds, as a list, a vector (named phase.id) numbering each major activity phase found in the record, a factor (named trip.act) labelling each row as being on-land, at-sea, or leisure at-sea activity. These two elements are as long as there are rows in tdr. This slot also contains two more vectors: one with the beginning time of each phase, and another with the ending time; both represented as POSIXct objects. See detPhase.

The slot dive.activity contains a data.frame, again with as many rows as those in tdr.consisting of three vectors named. dive.id, which is an integer vector, sequentially numbering each dive (rows that are not part of a dive are labelled 0), dive.activity is a factor which completes that in trip.act above, further identifying rows in the record belonging to a dive. The third vector in dive.activity is an integer vector sequentially numbering each postdive interval (all rows that belong to a dive are labelled 0). See detDive, and getDAct to access all or any one of these vectors.

dive, phases is a slot corresponding to a factor that labels each row in the record as belonging to a particular phase of a dive. See labDivePhase, and getDPhaseLab to access this slot.

land.thréshold, sea.threshold, dive.threshold, and speed.calib.coefs are each a single number representing parameters used for detecting phases, and calibrating the TDR. Except for the latter, these are mostly for internal use, and hence do not have an accessor function. See getSpeedCoef for accessing speed.calib.coefs.

The plot method for this class creates an interactive plot of the "TDR" object, labelling the identified observations with a chosen factor from the corresponding "TDRcalibrate" object. The argument surface allows for the inclusion/exclusion of identified surface readings.

TDRcalibrate-class 19

Objects from the Class

Objects can be created by calls of the form new ("TDRcalibrate", ...). The objects of this class contain information necessary to divide the record into sections (e.g. land/water), dive/surface, and different sections within dives. They also contain the parameters used to calibrate speed and criteria to divide the record into phases.

Slots

- tdr: Object of class "TDR", with concurrent time, depth, and possibly speed (if "TDRspeed").
 See Details.
- gross.activity: Object of class "list", must be the same as value returned by detPhase.
- dive.activity: Object of class "data.frame", must be the same as value returned by detDive.
- dive.phases: Object of class "factor", must be the same as value returned by labDivePhase.
- land.threshold: Object of class "numeric" the temporal criteria used for detecting periods on land that should be considered as at-sea.
- sea.threshold: Object of class "numeric" the temporal criteria used for detecting periods at-sea that should not be considered as foraging time.
- dive.threshold: Object of class "numeric" the criteria used for defining a dive.
- speed.calib.coefs: Object of class "numeric" the intercept and slope derived from the speed calibration procedure.

Methods

- attendance signature (obj = "TDRcalibrate", ignoreZ = "logical"): generates an attendance table for the TDR record; the duration of each dry and wet phase.
- getDAct signature(x = "TDRcalibrate", y = "missing"): extracts a data frame with vectors identifying all readings to a particular dive and postdive number, and a factor identifying all readings to a particular activity.
- getDAct signature(x = "TDRcalibrate", y = "character"): as the method for missing y, but selects a particular vector to extract.
- getDPhaseLab signature(x = "TDRcalibrate", diveNo = "missing"): extracts
 a factor identifying all readings to a particular dive phase.
- getDPhaseLab signature(x = "TDRcalibrate", diveNo = "numeric"): as the
 method for missing y, but selects data from a particular dive number to extract.
- extractDive signature(obj = "TDRcalibrate", diveNo = "numeric", id = "missing"):
 extract particular dives.
- getGAct signature(x = "TDRcalibrate", y = "missing"): extracts elements that
 divide the data into major wet and dry activities.
- getGAct signature(x = "TDRcalibrate", y = "character"): as the method for missing y, but extracts particular elements.
- $show \ \, \texttt{signature(object = "TDRcalibrate"): prints informative summary of the data.}$
- getTDR signature (x = "TDRcalibrate"); extracts the TDR object.

20 TDR-class

```
\begin{tabular}{ll} \textbf{getSpeedCoef} & \texttt{signature} (\texttt{x} = \texttt{"TDRcalibrate"}) \colon \textbf{extracts the speed calibration coefficients}. \end{tabular}
```

```
plot signature(x = "TDRcalibrate", y = "missing"): plot the TDR object, la-
belling identified sections of it (see Usage and Details).
```

Author(s)

```
Sebastian P. Luque (spluque@gmail.com)
```

See Also

TDR-class for links to other classes in the package

TDR-class

Classes "TDR" and "TDRspeed" for representing TDR information

Description

These classes store information gathered by time-depth recorders.

Usage

```
## S4 method for signature 'TDR, numeric, numeric':
extractDive(obj, diveNo, id)
```

Arguments

```
obj "TDR" object.
```

diveNo numeric vector or scalar with dive numbers to extract.

id numeric vector of dive numbers from where diveNo should be chosen.

Details

Since the data to store in objects of these clases usually come from a file, the easiest way to construct such objects is with the function readTDR for treireve all the necessary information. The methods listed above can thus be used to access all slots.

Objects from the Class

```
Objects can be created by calls of the form new("TDR", ...) and new("TDRspeed", ...).
```

TDR objects contain concurrent time and depth readings, as well as a string indicating the file the data originates from, and a number indicating the sampling interval for these data. TDRspeed extends TDR objects containing additional concurrent speed readings.

TDR-class 21

```
Slots
```

In class TDR:

```
file: Object of class "character", string indicating the file where the data comes from.
   dtime: Object of class "numeric", sampling interval in seconds.
   time: Object of class "POSIXct", time stamp for every reading.
   depth: Object of class "numeric", depth (m) readings.
   Class TDRspeed adds:
   speed: Object of class "numeric" speed (m/s) readings.
Methods
   as.data.frame signature (x="TDR"): Coerce object to data.frame.
   as.data.frame signature (x="TDRspeed"): Coerce object to data.frame.
   coerce signature (from="TDR", to="data.frame"): Coerce object to data.frame.
   coerce signature (from="TDRspeed", to="data, frame"); Coerce object to data, frame.
   getDepth signature (x = "TDR"): depth slot accessor.
   getDtime signature (x = "TDR"): sampling interval accessor.
   extractDive signature(obj = "TDR", diveNo = "numeric", id = "numeric"):
        extract particular dives.
   getFileName signature (x="TDR"): source file name accessor.
   plot signature (x = "TDR", y = "missing"): interactive graphical display of the data,
        with zooming and panning capabilities.
   plot signature (x = "TDRspeed", y = "missing"): As the TDR method, but also
        plots the speed slot.
   show signature (object = "TDR"): print an informative summary of the data.
   getTime signature (x = "TDR"): time slot accessor.
   getSpeed signature (x = "TDRspeed"): speed accessor for TDRspeed objects.
Author(s)
```

See Also

```
readTDR, TDRcalibrate-class.
```

Sebastian P. Luque (spluque@gmail.com)

22 zoc

zoc Interactive zero-offset correction of TDR data

Description

Correct zero-offset in TDR records, with the aid of a graphical user interface (GUI), allowing for dynamic selection of offset and multiple time windows to perform the adjustment.

Usage

Arguments

```
time
                  POSIXct object with date and time.
depth
                  numeric vector with depth in m.
                  known amount of meters for zero-offset correcting depth throughout the entire
offset
                  TDR record.
speed
                  numeric vector with speed in m/s.
                  vector of length 2, with lower and upper limits of time to be plotted.
xlim
phaseCol
                 factor dividing rows into sections.
xlab, ylab.depth, ylab.speed
                  strings to label the corresponding y-axes.
                 format string for formatting the x axis; see strptime.
xlab.format
```

Details

These functions are used primarily to correct, visually, drifts in the pressure transducer of TDR records. zoc calls plotDive, which plots depth and, optionally, speed vs. time with the possibility zooming in and out on time, changing maximum depths displayed, and panning through time. The option to zero-offset correct sections of the record gathers x and y coordinates for two points, obtained by clicking on the plot region. The first point clicked indicates the offset and beginning time of section to correct, and the second one indicates the ending time of the section to correct. Multiple sections of the record can be corrected in this manner, by panning through the time and repeating the procedure. In case there's overlap between zero offset corrected windows, the last one prevails.

Once the whole record has been zero-offset corrected, remaining points with depth values lower than zero, are turned into zeroes, as these are assumed to be values at the surface.

zoc 23

Value

zoc returns a numeric vector, as long as depth of zero-offset corrected depths.

plotDive returns a list with as many components as sections of the record that were zero-offset corrected, each consisting of two further lists with the same components as those returned by locator.

Author(s)

Sebastian P. Luque (spluque@gmail.com), with many ideas from CRAN package sfsmisc.

See Also

detDive

Index

```
*Topic arith
                                           .getSpeedCalib, 3
   diveStats, 10
                                           .getSpeedCalib
                                                  (diveMove-internal), 8
   rgPlot, 14
*Topic classes
                                           .getSpeedStats
   TDR-class, 20
                                                  (diveMove-internal), 8
   TDRcalibrate-class, 17
                                          as.data.frame.TDR-method
*Topic datasets
                                                  (TDR-class), 20
   sealMK8.15
                                          attendance. 9
*Topic hplot
                                          attendance (TDRcalibrate-class),
   raPlot. 14
*Topic internal
                                          attendance, TDR calibrate, logical-method
   diveMove-internal.8
                                                  (TDRcalibrate-class), 17
*Topic iplot
                                          austFilter,1
   zoc, 22
*Topic iteration
                                          calibrateDepth, 3,9
   austFilter, 1
                                          calibrateSpeed, 9
*Topic manip
                                          calibrateSpeed(calibrateDepth), 3
   austFilter.1
                                          coerce, TDR, data.frame-method
   calibrateDepth.3
                                                  (TDR-class), 20
   detDive.4
                                          createTDR (readTDR), 13
   detPhase, 5
   distSpeed, 7
                                          data.frame, 11
                                          detDive, 3, 4, 4, 6, 18. 19. 23
   doSpeedCalib, 16
                                          detPhase, 3, 4, 5, 5, 12, 19
   readLocs, 12
                                          distSpeed, 2, 7
   readTDR, 13
                                          diveMove (diveMove-package), 9
   raPlot. 14
                                          diveMove-internal, 8
*Topic math
                                          diveMove-package, 9
   calibrateDepth, 3
                                          diveStats, 10, 15
   distSpeed, 7
                                          doSpeedCalib, 3, 4, 15, 16
   diveStats, 10
   doSpeedCalib, 16
                                          extractDive(TDRcalibrate-class),
*Topic package
   diveMove-package, 9
                                          extractDive, TDR, numeric, numeric-method
.cutDive(diveMove-internal), 8
                                                  (TDR-class), 20
.descAsc(diveMove-internal).8
                                          extractDive, TDRcalibrate, numeric, missing-method
.diveIndices (diveMove-internal),
                                                  (TDRcalibrate-class), 17
.getInterval(diveMove-internal),
                                          getAct (detPhase), 5
                                          getDAct. 18
```

INDEX 25

getDAct(TDRcalibrate-class), 17	locator, 23
getDAct, TDRcalibrate, character-method	
(TDRcalibrate-class), 17	matrix, 16
getDAct, TDRcalibrate, missing-method	
(TDRcalibrate-class), 17	plot, TDR, missing-method
getDepth, 18	(TDR-class), 20
getDepth (TDR-class), 20	plot, TDRcalibrate, missing-method
getDepth, TDR-method (TDR-class),	(TDRcalibrate-class), 17
20	plot, TDRspeed, missing-method
getDive(diveStats), 10	(TDR-class), 20
getDPhaseLab, 18	plotDive, 18
getDPhaseLab	plotDive(zoc), 22
(TDRcalibrate-class), 17	postscript, 16
(IDRCaliblate-Class), I/	had
getDPhaseLab, TDRcalibrate, missing-met	readLocs, 12
(TDRcalibrate-class), 17 getDPhaseLab, TDRcalibrate, numeric-met	readTDR, 13, 15, 18, 21
getDPhaseLab, TDRcalibrate, numeric-met	rmsDistFilter(austFilter), 1
(TDRCalibrate-class), 17	rqPlot, 14, 17
getDtime(TDR-class), 20	-4-13,13,0
getDtime, TDR-method (TDR-class),	sealMK8,15
20	show, TDR-method (TDR-class), 20
getFileName (TDR-class), 20	show, TDRcalibrate-method
getFileName, TDR-method	(TDRcalibrate-class), 17
(TDR-class), 20	stampDive.9
getGAct (TDRcalibrate-class), 17	stampDive(diveStats), 10
getGAct, TDRcalibrate, character-method	strptime, 8, 12, 13, 22
(TDRcalibrate-class), 17	strptime, 0, 12, 13, 22
getGAct, TDRcalibrate, missing-method	TDR, 18, 20
(TDRcalibrate-class), 17	TDR-class, 9, 20
getSpeed, 18	
getSpeed (TDR-class), 20	TDR-class, 20
getSpeed(TDR C1333), 20 getSpeed, TDRspeed-method	TDRcalibrate, 18
(TDR-class), 20	TDRcalibrate-class, 3, 10, 12, 17, 21
	TDRcalibrate-class, 17
getSpeedCoef, 18	TDRspeed-class (TDR-class), 20
getSpeedCoef	track(distSpeed),7
(TDRcalibrate-class), 17	
getSpeedCoef, TDRcalibrate-method	zoc, 3–5, 12, 22
(TDRcalibrate-class), 17	
getTDR, 18	
getTDR(TDRcalibrate-class), 17	
getTDR, TDRcalibrate-method	
(TDRcalibrate-class), 17	
getTime, 18	
getTime(TDR-class), 20	
getTime, TDR-method (TDR-class), 20	
grpSpeedFilter(austFilter), 1	
2 T - T	
labDive (detDive), 4	
labDivePhase, 18, 19	
labDivePhase (detDive), 4	